



DID SABZ



**Dispenser Plus**

**DS 4LM+**

# Dispenser Plus

DS 4LM +




شرکت دید سبز

یک سال گارانتی / ده سال پشتیبانی فنی

## Technical Data:

<b>System Control :</b>	<b>Microcontroller</b>
<b>Dimensions (W x D x H):</b>	300 x 520 x 290 mm
<b>Weight:</b>	approx. 18 kg
<b>Voltage:</b>	220 v 50 Hz
<b>Paraffin Capacity:</b>	Up to 2.3 Liters

- Holding three warm environments:  
A Large Tank for melting paraffin, another Tank for heating the Plastic Molds and a Hot Plate
- Online temperature display with 0.1°C accuracy for tank and Hot plate
- Equipped with an electric valve which controls the flow volume
- Delay timer which could be set up to 10 days for holidays
- The alarm lights which conform to each function of the device
- The electric valve of the timer of which could be set from 0 to 300 seconds and also could be turned off
- Two totally separate temperature controlling environments
- Controlling the paraffin tanks using the thermal sensors from 30°C to 80°C
- Controlling the Hot Plate temperature from 30°C to 80°C using the thermal sensors
- Ease of reading the temperature via the flat panel and LCD display

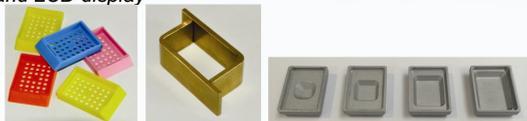

In Conformity to IEC61010-1

## مشخصات و ویژگی ها:

- سیستم کنترل دستگاه از نوع میکروکنترلر
- راحتی قرائت دمای پانل فلت با نمایشگر LCD
- مجهز به پدال پایی برای آزاد بودن هر دو دست اپراتور جهت بالا بردن دقت و سرعت عمل قالب گیری
- سه محیط کنترل دما بصورت مجزا:
- یک مخزن بزرگ استیل برای ذوب پارافین، یک مخزن کوچک استیل برای گرم کردن تیشوتک و یک صفحه گرم
- کنترل قطع و وصل سیال توسط شیر الکتریکی بدون چکه و با قطعی کامل
- مجهز به چراغ زیر شیر الکتریکی
- جلوگیری از انسداد مسیر خروجی پارافین با تامین حرارت لازم
- کنترل مقدار حجم خروجی سیال با استفاده از یک سیستم الکترونیکی
- ظرفیت مخزن بزرگ ۲/۳ لیتر، ظرفیت مخزن کوچک ۰/۶ لیتر و ابعاد صفحه گرم ۲۳ \* ۲۵ سانتی متر
- تایمر تاخیری قابل برنامه ریزی برای ایام تعطیلات تا ۱۰ روز
- نصب چراغ های هشدار دهنده متناسب عملکرد دستگاه
- شیر برقی قابل کنترل توسط تایمر صفر تا ۳۰۰ ثانیه که می توان تایمر را خاموش کرد.
- کنترل صفحه گرم در محدوده دمایی از ۳۰ الی ۸۰ درجه سانتیگراد با استفاده از سنسور حرارتی

DID SABZ Co.  
ADD.: Azadi Blvd., Opposite To Customs House  
Ostad Tanha Ave., Unit 6, Urmia-Iran  
Postal Code: 57176-66593  
Tel-Fax: (+98) 44 32771214-32781173



www.didsabz-co.com

E-mail: info@didsabz-co.com

شرکت دید سبز (سهامی خاص)

آدرس: ارومیه، بلوار آزادی، رویروی گمرک،  
خیابان استاد تنها، واحد ششم، دست راست  
شرکت دید سبز کد پستی: ۵۷۱۷۶-۶۶۵۹۳  
تلفن و دورنویس: ۰۴۴) ۳۲۷۷۱۲۱۴-۳۲۷۸۱۱۷۳